

BODY-JET RESULTS

Patients will demand the minimally invasive Body-Jet due to the smooth result, quick recovery and minimal discomfort. Body-Jet can be used under local anesthesia, further minimizing risk for the patient.



“The innovative Body-Jet affords physicians a multitude of benefits and is creating exciting opportunities for fat grafting. My practice is benefitting greatly from increased booking, referrals and satisfied Body-Jet patients.”
— Barry DiBernardo, MD



“The Body-Jet has opened the doors to greater possibilities for routine fat and stem cell harvesting and transfer while vastly improving everyday liposuction results. Our creativity and artisanship are now no longer limited by our technology.”
— Todd Malan, MD



Courtesy of Barry DiBernardo, MD



Courtesy of Todd Malan, MD

Courtesy of Jeffrey Caruth, MD

Experience the difference today. Realize the rewards of the natural flow of water for aesthetic body contouring. Visit us at www.eclipsemed.com to learn more and be part of this new cutting edge medical technology.



body-jet®
the gentle way of liposuction

The Body-Jet system includes:

- ✦ WAL Applicator, Sterile, Single-use
- ✦ Infiltration Cannulas (3)
- ✦ Irrigation and Aspiration Cannulas (14)
- ✦ Suction bags (24)
- ✦ Cannula cleansing brush (2)

Body-Jet specifications

DIMENSIONS (cm) 178H x 542W x 60D
WEIGHT 101lbs (46kg)

The AquaShape®FT LipoCollector II System includes:

- ✦ LipoCollector Container
- ✦ LipoCollector Disposable Set (1 set)
- ✦ Mesh Filter, 200µm, sterile (1)
- ✦ Mesh Filter, 250µm, sterile (1)
- ✦ Mesh Filter, 315µm, sterile (1)
- ✦ Connection tubing to suction container, 3m, sterile (1)

AquaShapeFT LipoCollector NOTE: For optimal harvesting of fat during a liposuction, a slight modification of the surgical technique (speed, cannula size and suction force) may be needed. Testing of the device in a 'liposuction' only case first is highly recommended. Eclipsemed, Ltd. and Human Med AG do not take any responsibility for the quality of the harvested material or surgical procedure outcome.

Body-Jet® is a Class II product. The Body-Jet® has patents in the United States (10/520, 209; 10/523, 294; 10/322, 524), other patents and patents pending. Body-Jet is a registered trademark of Human Med AG, Schwerin, Germany.

CAUTION: Federal law restricts sale of Body-Jet by or on the order of a physician.

Manufactured by



Distributed in the United States by Eclipsemed, Ltd.

16850 Dallas Parkway, Dallas, TX 75248

800-759-6876 • www.eclipsemed.com • www.mybodyjet.com



Water-jet assisted
lipoplasty for Aesthetic
Body Contouring

Transform your practice
and patients with the
natural flow of water




THE REVOLUTION STARTS NOW.

As the industry's only
water-jet assisted
liposuction, Body-Jet is
setting a new standard
and transforming the
way doctors and patients
view liposuction. Using a
gentle, pulsating spray
of fluid, Body-Jet offers
a new alternative to
conventional techniques.

body-jet[®]
the gentle way of liposuction


BODY-JET

Water-jet assisted liposuction (WAL) with Body-Jet activates a slightly pressurized fan-shaped laminar jet to anesthetize tissues. Then a customized cannula loosens and removes the fat simultaneously so that functional tissue structures are spared. Fat is then harvested and collected in the AquaShape[®]FT LipoCollector II.

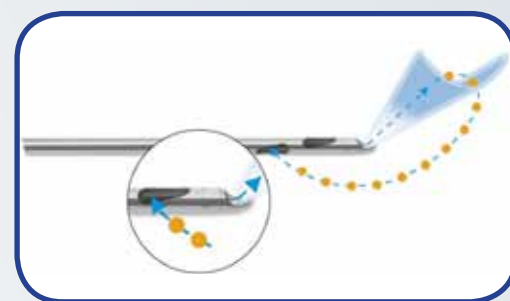

This unique process involves two phases:

- *Infiltration of analgesics and adrenaline* — usually 50% - 70% less TLA infiltration volume and exposure time versus conventional methods
- *Simultaneous irrigation and aspiration* — fat cells are easily 'released' by the flow of fluid through the cannula irrigation

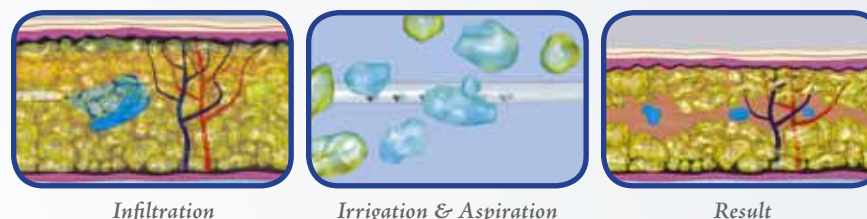

TECHNOLOGY TO TRANSFORM YOUR PRACTICE

For physicians, the benefits speak for themselves and ultimately result in increased revenues and satisfied patients.

SAFETY

- Minimal procedure risk
- Less intra-operative tissue trauma resulting in less swelling
- Reduced exposure to tumescent fluid

TIME

- Decreased waiting time
- Shorter overall procedure time
- Minimal surgeon effort due to gentle method

RESULTS

- Accurate contours and better intra-operative assessment
- Reduced correction rate
- Quick and efficient fat harvesting technique

AQUASHAPE



AQUASHAPE[®]FT LIPOCOLLECTOR II

The AquaShape[®]FT LipoCollector II is suitable for the collection of fat intended for re-injection (fat grafting/fat harvesting) or laboratory purposes. The aspirate collected from a Body-Jet liposuction procedure usually contains very little blood. The LipoCollector separates fat and fluid instantly after aspiration, leaving only a very short time for the fluid to affect the lipocytes. The container is fully autoclavable and does not have to be opened.

BENEFITS

- Collect up to one liter of fat during a liposuction procedure
- Gentle handling of fat cells
- No need for centrifuging
- Fully autoclavable device

