

tixel[®]

Peer Reviewed, Proven Success
Fractional Skin Rejuvenation





Tixel is a fractional skin rejuvenation system designed to treat stubborn skin conditions (lip lines and other wrinkles) as well as the delicate skin on your eyelids, neck and hands.

Powered by Thermo-Mechanical Action® (precisely controlled heat), incredible results can be achieved by simply applying a unique patented titanium tip technology to targeted areas.

Skin emerges smoother, tighter and younger-looking.

This sophisticated technology is suitable for all skin types and all ages and can be safely used year-round in any climate.





Treatment Areas

Forehead lines (worry lines)

Frown lines (glabellar lines)

Eyelid wrinkles

Crow's feet (periorbital lines)

Swollen mounds under the lower eyelid (festoons)

Smile lines (nasolabial folds)

Lip lines

Creased chin

Neck lines

Arms

Décolleté (chest)


Stretch marks

Hands

and more




How Tixel Works

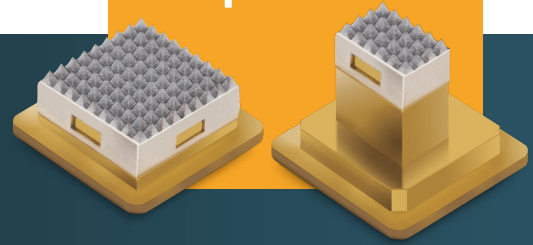


The titanium tip is heated to a precise temperature with thermal energy. Then the handpiece is moved across the treatment areas, briefly contacting the skin and providing uniform coverage. The extent of energy transfer is determined by the tip's exposure time to the skin. The longer the pulse duration, the greater the amount of energy transferred. This causes a cascade of rejuvenation: new collagen, new elastin, fibroblasts, and more.


Pulse Duration

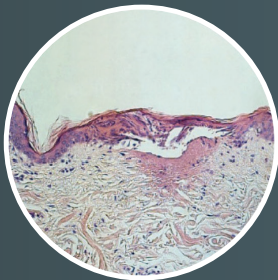

SHORT 6-8 ms

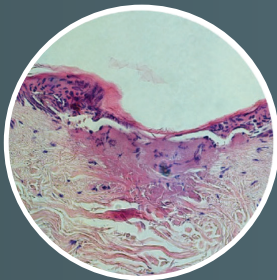

MEDIUM 10-12-14 ms

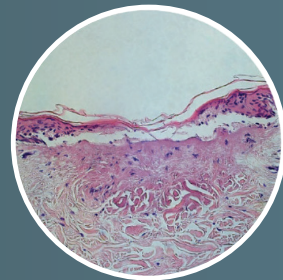

LONG 16-18 ms

After treatment, patients can expect:

- ✓ Slightly swollen and pink skin for 3 to 4 hours
- ✓ Mild erythema and heat sensation for 2 to 3 hours
- ✓ Downtime from 0 to 2 days, depending on treatment intensity
- ✓ Micro-crusting disappears in about 3 days

The simplicity and elegance of Tixel is undeniable.

Simplicity

- ✓ No anesthesia; analgesic creams are optional
- ✓ Low pain
- ✓ Non-invasive, no bleeding
- ✓ Minimal downtime (0-2 days for most patients)
- ✓ No eyewear needed for patients or clinician

Elegance

- ✓ Treats skin conditions from head to toe
- ✓ Suitable for all skin types and all ages
- ✓ Can be used year-round, even in summer
- ✓ Radiation-free, safe and easy for patients and clinicians


Before & After Results



Number of Tixel treatments 5



Number of Tixel treatments 3

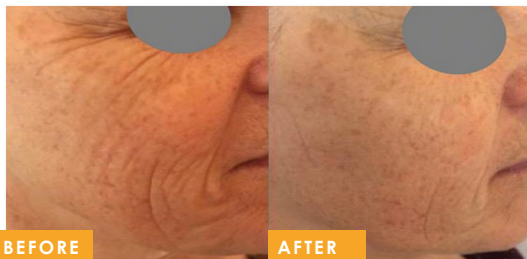

Number of Tixel treatments 3
Tixel settings: 8-400, 10-500, 12-400



Number of Tixel treatments 3
Tixel settings: 14-700




FIND YOUR PERFECT FIT.

See how the classic Tixel C stacks up against the next-generation Tixel 2
—and discover which one belongs in your treatment room.

FEATURES



Simultaneous pre-warmed & ready	2	1
Repeat mode	Yes	N/A
Hot swap of handpieces	Yes	N/A
Compatible with large handpiece	Yes	Yes
Compatible with precision handpiece	Yes	Yes
Start up time	-3 minutes	6 minutes
Cool-down (shut-down)	-3 minutes	6 minutes
Operator notifications	Full screen pop-up window	Flashing icon
Smaller precision handpiece included for safe, targeted treatment of delicate areas such as the eyes and mouth — plus a second handpiece for larger treatment zones	Yes	N/A (small handpiece available for purchase separately)

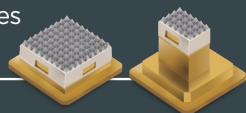

Users should always consult the Operator Manual for recommended settings.

Tixel, Tixel 2, Thermo-Mechanical Action, and Novoxel, Simple Elegance are all registered trademarks of Novoxel Ltd.

Copyright 2024, Novoxel Inc. All rights reserved.

Contact us: www.clarionmedical.com info@clarionmedical.com @Clarion_medical